

ALE TOTOTAL

(3)

DTIC FILE COP'

TECHNICAL REPORT RD-PR-86-3

SELECTIVE DETECTION OF MINIMUM SMOKE STABILIZERS: FINAL REPORT

STIC ELECTE MAY 2 8 1987

Dr. James G. Carver
Propulsion Directorate
Research, Development and Engineering Center

JULY 1986



U.S. ARMY MISSILE COMMAND

Redstone Arsenal, Alabama 3500-50

Cleared for public release; distribution unlimited.

DISPOSITION INSTRUCTIONS

DESTROY THIS REPORT WHEN IT IS NO LONGER NEEDED. DO NOT RETURN IT TO THE ORIGINATOR.

DISCLAIMER

THE FINDINGS IN THIS REPORT ARE NOT TO BE CONSTRUED AS AN OFFICIAL DEPARTMENT OF THE ARMY POSITION UNLESS SO DESIGNATED BY OTHER AUTHORIZED DOCUMENTS.


TRADE NAMES

USE OF TRADE NAMES OR MANUFACTURERS IN THIS REPORT DOES NOT CONSTITUTE AN OFFICIAL INDORSEMENT OR AFPROVAL OF THE USE OF SUCH COMMERCIAL HARDWARE OR SOFTWARE.

REPORT DOCUMENTATION PAGE OMB NI Exp. Da		Form Approved OMB No 0704 0188 Exp. Date: Jun 30: 1986			
14 REPORT SECURITY CLASS FICATION	16 RESTRICTIVE MARKINGS				
Unclassified 3 DISTRIBUTION AUTHORITY 3 DISTRIBUTION AVAILABILITY OF RE		JE DEDOGT			
24 SECURITY CLASSIFICATION AUTHORITY		3 DISTRIBUTION	AANITARITITA (JE MEPUKI	
26 DECLASSIFICATION DOWNGRADING SCHEDULE					distribution
4 PERFORMING ORGANIZATION REPORT NUMBE	R(S)	5 MONITORING	ORGANIZATION	REPORT NU	MBER(S)
RD-PR-86-3	9-PR-80-3				
64 NAME OF PERFORMING ORGANIZATION	66 OFFICE SYMBOL	7a NAME OF M	ONITORING ORG	ANIZATION	
Commander	(If applicable)				
L. S. Army Missile Command 6c ADDRESS (City State and 2/P Code)	AMSMI-RD-PR	76 ADDRESS (Ci	ty State and 716	Code)	
or mount is saily, state and air code)] · · · · · · · · · · · · · · · · · · ·	iy, Jiait, difu Lif	-WE	
Redstone Arsenal, AL 35898-5	249				
Ba NAME OF FUNDING SPONSORING	86 OFFICE SYMBOL	9 PROCUREMEN	T INSTRUMENT I	DENTIFICAT	ION NUMBER
ORGANIZAT ON	(If applicable)				
Bc ADDRESS (City State and 21P Code)	<u> </u>	10 SOURCE OF	UNDING NUMBE	RS	
		PROGRAM	PROJECT	TASK	WORK UNIT
		ELEMENT NO 72KD01	00	NO O	ACCESSION NO 7280
(Include Security Classification)		/2KD01	1 00	1	7280
Selective Detection of Minimum	m Smoke Stabilia	zers: Final	Report		
12 PERSONAL A THORIS					
Carlet, Limes C.					
	y 8+ 10 Sep 85	14 DATE OF REPO July 1986			PAGE COLS.
The speed of the state of the s	ject has been ac				
nine, reddit, er frett Lestine Lemante elletin	vre equipment (i	in mechanica	l, chemical	. or no	ndestructive
Tradition to mante controlled	TA SHARE THE TOWN				
	tabilizers,	. Minimum Se	ке, Double	Base,	Propellest,
├ +	what sais. Il				
TO SERVICE TO THE PROPERTY.	in a dent ty by brick o	umber			
	to the contract of				
the provide the state of the st	$\label{eq:continuous} \begin{array}{cccccccccccccccccccccccccccccccccccc$	1 1 1	in district the second second	ed Land	in the Baresmith
		€ 1 × 31			
		I			
			7.		
t e e			1		u. C

TABLE OF CONTENTS

		Page
LIST	OF FIGURES	iv
ι.	INTRODUCTION	t
II.	EXPERIMENTAL	2
III.	RESULTS	3
ī۷.	CONCLUSIONS	12

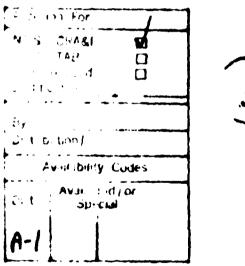

RE: Report has Military Application, Distribution Statement A Distribution Statement A is Correct for this report.

Per Dr. James G. Carver, Army Missile, AMSMI-RD-PR

I. INTRODUCTION

The Armed services lately have been expressing an increased interest in the safe life of propellants. The Army has an Independent Review Team examining the status of gun and missile propellants, in storage, for present safety and remaining safe life. The Army has also placed a high priority on the evaluation of the safety of the M55 chemical rocket. There are over 600,000 rounds stored in several locations around the U.S. and in the Pacific which are planned for disposal. However, before these rounds can be shipped to a central location for destruction, it is necessary that their safety be evaluated.

Of the 21 Army missile systems currently in development or fielded, 12 use double base, composite modified double base, or minimum signature propellants as the main grain. Most of the 21 systems use 1 of these propellants in their ignitors or gas generators. Many other systems with these propellants are in storage under the demil 11-P account awaiting disposition.

The safe life as well as the usable service life of these missiles is directly dependent upon the levels of the various stabiliziers and their degradation products present in the propellant. Therefore, it is critical that these levels he accurately known at all times. Gas chromatography (GC) and liquid chromatography (LC) are the two procedures currently used in the propellant industry to determine stabilizer levels. The most frequently used detector in GC is the Flame ionization Detector (FID), and the most favored LC detector is the LV detector. Both of these systems are used because of their ability to detect almost overse ingredient in the propellant. This is a problem, however, with the complex mixtures being used in the newer minimum signative propellant formulations. It can be difficult to completely resolve the stabilizers from the other ingredients in these complex mixtures.

A liquid chromatorrophic detector, that has received some attention in forensic science, is the electrochemical detector. This detector has the ability to be selective and very sensitive. In the oxidative mode, it can be adjusted to detect only aromatic amines and phenois at ng levels. The most trequently used stabilizers in the United States are of this type. Therefore, a study to evaluate the use of the electrochemical detector in the detection of propellant stabilizers was conducted.

II. EXPERIMENTAL

Experiments were conducted using the liquid chromatographic system which consisted of a Waters Associates 720 Controller, an M6000A pump, and a Wisp 710B Automatic Sample Injector, a Perkin-Elmer LC85 variable wavelength UV detector with autocontrol, and a BAS LC-4B dual amperometric detector.

Separations were performed on a DuPont Zorbax 5 µm ODS 25 cm X 4.6 mm column or a Waters C-18 NOVAPAC 5 µm column. In both cases, a Brownlee Labs 3 cm X 4.6 mm ODS Guard Column was placed in front of the analytical column. The eluting solvents compositions varied with the propellant, but were mixtures of acetonitrile, methanol, tetrahydrofuran, and 0.1 M sodium acetate adjusted to a pH of 5 with perchloric acid. All ingredients were reagent grade. Solvents were UV grade HPLC quality, and water was distilled, deionized, and filtered.

The UV detector was typically set to 214 nm. The electrochemical detector could be arranged so that the working electrodes were parallel or in series. Most of this work was performed with the series conformation. The working electrodes were dual glassy carbon. The potentials were set relative to a silver/silver chloride reference electrode.

The analog output from each detector was digitized and sent to a Hewlett-Packard 3357 lab automation system for integration and data reduction.

Samples of propellant were prepared in two similar procedures. The propellant was hand chopped to a coarse powder. In one procedure the propellant powder (25 to 100 mg) is accurately weighed into a 10 mL test tube. The internal standard is then accurately weighed into a small aluminum pan and placed in the test tube. Five to seven mL of a 50 percent ACN solution is then added and the sample is sonicated with a microprobe for 10-15 min. In the alternate procedure 25 mg of the propellant powder was accurately weighed into a 5 mL vial. The internal standard was accurately weighed into a small aluminum pan and placed in the vial. Four milliliters of dichloroethane was pipeted into the vial. The vial was capped and the propellant soaked overnight.

Once the propellant had been extracted by one of the above procedures, the supernatent liquid was pipeted into a 5 mL syringe fitted with a 0.45 .m disposable filter and filtered into a sample vial.

ዸቔዿኯቔ፟፟ኯቔኇቔቔዿቑዸቔዸቔቔቔቔቔቔ፟ዸኯዾዾዄቔዾዾዸቜቔዄጚኯ፟ዄዄጜኇዿኇዾዹዄኯዄኇቜኇፚፙዄዄዄዀጜኯዾ፞ፙቜዄዄጜዀዄኯዹ፟ኯዄቔዹቔቔቔ

III. RESULTS

One of the problems encountered when selecting the conditions to be used with an electrochemical (EC) detector was the choice of the mobile phase. For the detector to function properly, the eluting solvent must conduct electricity. This usually requires a minimum ion concentration of 0.05 M in the solvent. Thus, for a 50 percent organic phase solution the aqueous phase must be 0.1 M. The ions used are typically buffered salts of acetic, phosphoric, or citric acids. These give pH values within the safe operating limits of silica based columns. If a solution is not consumed in a few days, it should be discarded since it can act as a growing medium for biological species. Sodium acetate was selected for this study since it has a low UV cutoff and solutions have a longer shelf life than the phosphate and citrate salts.

Solutions prepared from sodium acetate and glacial acetic acid were found to give noisy background and a high UV cutoff. A solution of 0.1 M sodium acetate adjusted with perchloric acid to a pH of 5 was found to give the best results.

Solutions of the 32 propellant ingredients listed in Table 1 were prepared in acetonitrile. The EC detector was set to its highest recommended exidative potential of $1.20~\rm V$ and each ingredient was eluted through the system. As shown in Table 1, only the stabilizers which have an amine or hydroxyl on an aromatic ring are detected in the exidative mode at $1.20~\rm V$. As a point of interest, when the upstream electrode is set to $-0.9~\rm V$, (the reductive mode) and the down stream electrode to $\pm 1.0~\rm V$ most of the materials with nitro groups could also be detected.

The response of the EC detector to a substance is dependent on the concentration of the substance, the half-wave potential of the substance and the potential of the working electrode. For a given concentration of a substance the response vs potential curve looks much like a pH curve for titration of a monoprotic acid. There is a portion of the curve below the threshold potential where no response above the background current is detected. Above the threshold potential the response increases as the oxidative potential is increased until a plateau is obtained. Applying higher potentials generally has little effect on the response to the ingredient but does increase the background current. The potential midway up the response curve is the halfwave potential.

The ingredients which responded to the EC detector at 1.20 V were analyzed at a series of oxidative potentials between 0.6 V and 1.20 V. The normalized response vs potential curves are given in Figures 1, 2, and 3. As shown in these figures, with careful selection of the oxidative potential it is possible to detect the parent stabilizer and not its degradation products. For example, in Figure 3 at 0.80 V diphenylamine is detected while nitroso DPA, 2-NDPA, 4-NDPA, and the dimitro and trinitro derivatives do not respond in the electrochemical detector. It is also evident that if the original stabilizer is 2-NDPA it can be detected at 0.9 V or 1.0 V while the dimitro-and trinitro-DPA derivatives do not respond. This might become useful if resolution becomes a problem.

TABLE 1. Detectability of Propellant Ingredients

Ingredient	UV (214 nm)	EC (1.20 V)
Resorcinol	Yes	Yes
2 Nitro Resorcinol	Yes	Yes
2,4 dinitro Resorcinol	Yes	Yes
Triacetin	Weakly	No
HMX	Yes	No
RDX	Yes	No
DEGDN	Yes	oli
Nitroglycerin	Yes	No
BTTN	Yes	No
MNA	Yes	Yes
Nitroso MNA	Yes	Yes
Acaradite I	Yes	Yes
Acaradite II	Yes	Yes
2,4 DNT	Yes	No
2,6 DNT	Yes	No
2,4,6 TNT	Yes	No
Tetryl	Yes	No
TEGDN	Yes	No
TMETN	Yes	No
Methyl Centralite	Yes	Yes
Ethyl Centralite	Yes	Yes
Diphenylamine	Yes	Yes
Nitroso-DPA	Yes	Yes
2-NDPA	Yes	Yes
4-NDPA	Yes	Yes
2,2' DN-DPA	Yes	Yes
2,4' DN-DPA	Yes	Yes
2,4,4* TH-DPA	Yes	Yes
Methyl phthalate	Yes	No
Lthy! whthalate	Yes	No
Propyl phthalate	Yes	No
Propyl Adipate	Weakly	No

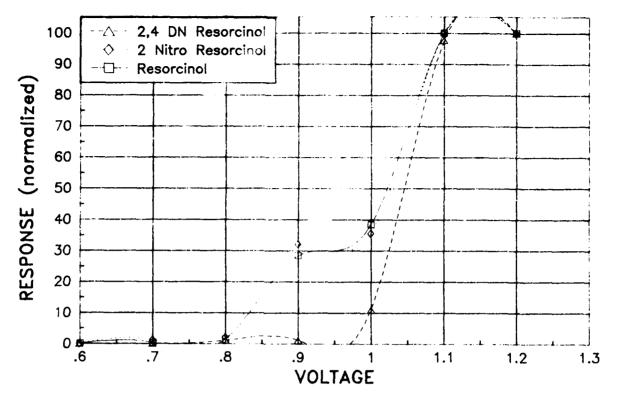


Figure 1. EC data - resorcinol.

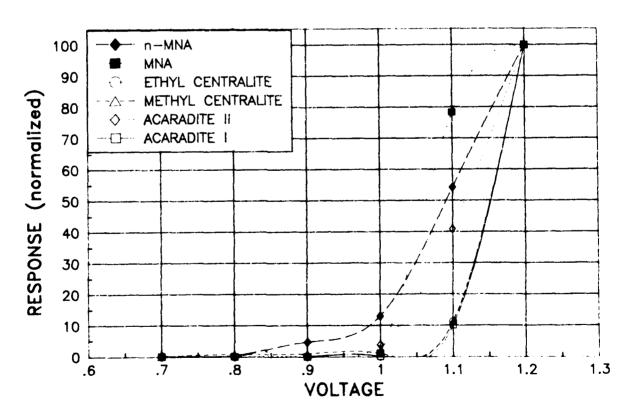


Figure 2. EC data - miscellaneous.

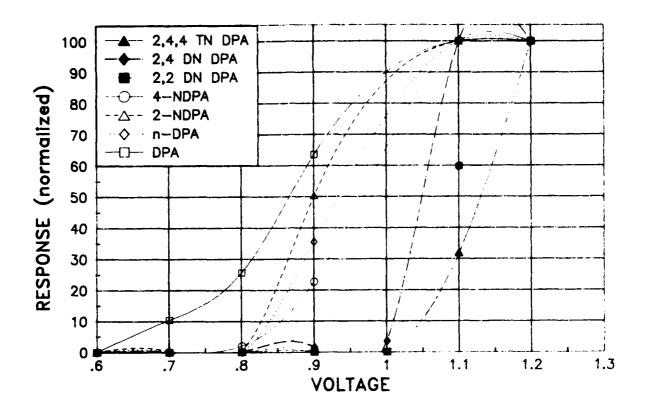


Figure 3. EC data - DPA.

The sensitivity of the EC detector to the stabilizers and their products was examined. Solutions of the 16 materials were prepared at 5 concentrations 1 mg/mL, 100 µg/mL, 10 µg/mL, 1 µg/mL and 0.1 µg/mL. The last concentration represents 0.001 percent in the propellant sample prepared by standard procedure. The results are shown in Figures 4, 5, and 6. The lowest detectable amount of material was assumed to have a response or peak height of twice the background noise level. Table 2 indicates the estimated minimum detectable amount of each material. The response of the detector is linear up to about 750 nA on the highest setting. Table 2 also indicates the upper limit of linearity for each compound. During this study, good results were not obtained from 2,4 Dinitro resorcinol. After the experiment was completed, it was found that a small baseline disturbance occurred at the same time that the material eluted. This disturbance was probably an injection phenomenon.

Comparison of the lower limits of detectability for the UV detector at 214 nm with the lower limits for the EC detector indicated that the EC was much more sensitive than the UV at 214 nm. This can be a problem if the sample is too concentrated. Originally, during this analyses poor accuracy was obtained because the linearity of the EC detector was exceeded. By preparing samples with one tenth to one half the concentrations typical for UV detection, this problem can be avoided.

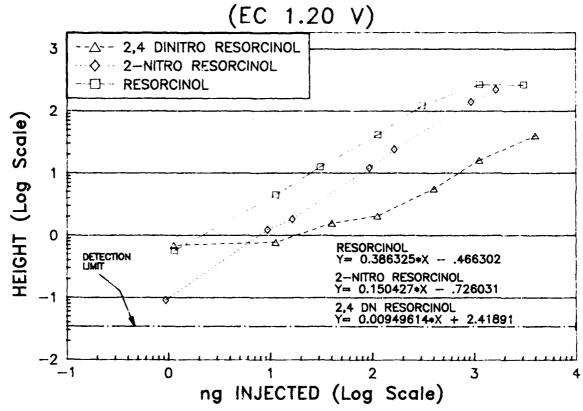


Figure 4. Propellant sample results - resorcinol.

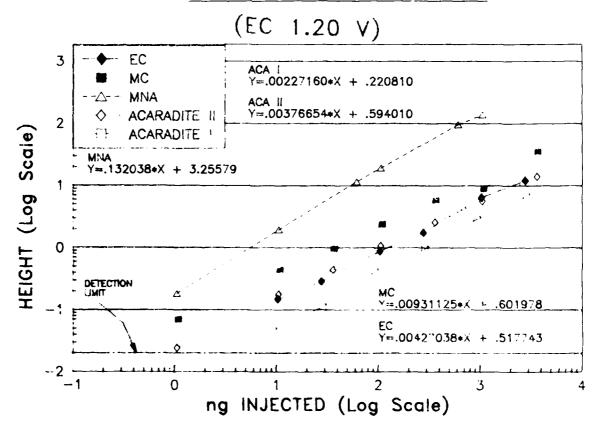


Figure 5. Propellant sample results - miscellareous.

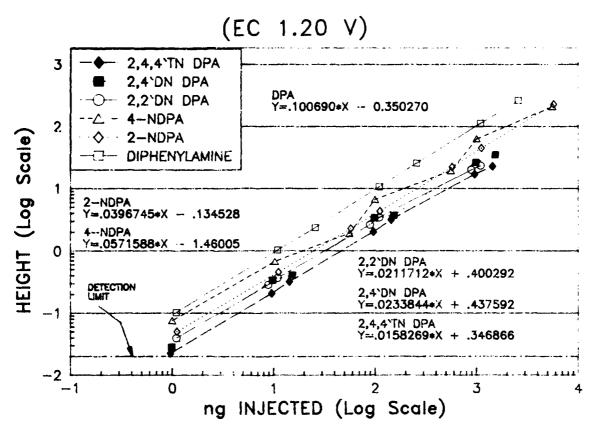


Figure 6. Propellant sample results - DPA.

TABLE 2. Linearity and Detection Limits

Material	Linearity (EC 1.2 V)	Detection (EC 1.2 V)	Detection (UV 214 nm)
Resorcinol	260 ng	.05 ng	10 ng
2-Nitro Resorcinol	670 ng	•4 ng	80 ng
Acaradite I	44000 ng	5.0 ng	15 ng
Acaradite II	26000 ng	1.0 ng	20 ng
MNA	730 ng	·1 ng	25 ng
MC	10000 ng	•3 ng	35 ng
EC	23000 ng	1.0 ng	50 ng
DPA	990 ng	•2 ng	20 ng
2-NDPA	2500 ng	•5 ng	32 ng
4-NDPA	1800 ng	.3 ng	58 ng
2,2'DNDPA	4700 ng	•6 ng	10 ng
2,4'DNDPA	4300 ng	.8 ng	50 ng
2,4,4'TNDPA	6300 ng	1.0 ng	15 ng

During the course of this study, it was found that the precision and accuracy of propellant analyses were poor when compared to GC or LC-UV techniques. A study of the sources of error disclosed the following:

- Use of methanol in the eluting solvent, even at less than 10 percent concentrations, slowly desensitizes the glassy carbon electrodes. This changes the response of the detector to the materials. The effect can cause significant changes in calculated stabilizer levels over a period of even 1 hour.
- Flow reproducibility is very critical since the EC detector is a concentration sensitive detector like UV rather than a mass sensitive detector like FID. Higher flow rates, such as 2 to 4 mL/min, were found to result in greater precision than slow flow rates such as 0.5 to 1 mL/min. Proper solvent preparation such as, filtering, mixing, and degassing are also critical to reproducible retention time and area precision.
- When real world samples were used it was found that the propellant matrix can contain ingredients that slowly desensitize the glassy carbon electrodes. The rate is dependent on the sample concentration, volume injected, and the number of injections made. Typically, this required the electrodes to be cleaned once a week, a 5-minute operation. Heavier use, as would be expected in an active QA lab, might require electrode cleaning every other day.
- The precision and accuracy of the procedure was evaluated by analyzing a sample of HYDRA 70 propellant. Figures 7 and 8 represent the typical responses of the UV detector at 214 nm and EC detector at 1.05 V in the chromatographic analysis of HYDRA 70 propellant. Table 3 indicates the results of triplicate analyses of two samples of RAD84K003-018. Within sample precision is good with a standard deviation of 0.019 percent and 0.029 percent and a relative standard deviation of 0.98 percent and 1.44 percent. Between sample precision is also good with a standard deviation of 0.023 percent and a RSD of 1.14 percent. The accuracy is also excellent. Gas chromatographic analysis found a 2-NDPA level of 1.98 percent + 0.01, while the LC/EC found 2.00 percent + 0.02.

Capy available to Diffe dose and permit fully legible reproduction

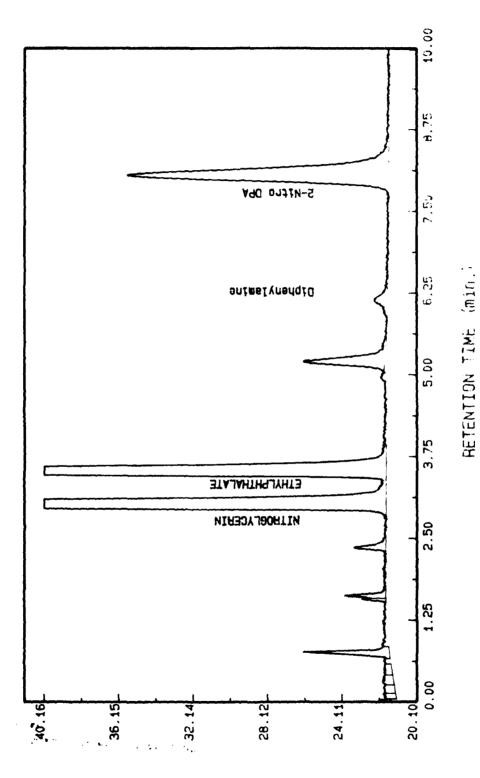


Figure 7. Typical responses of LC/UV detector at 214 nM.

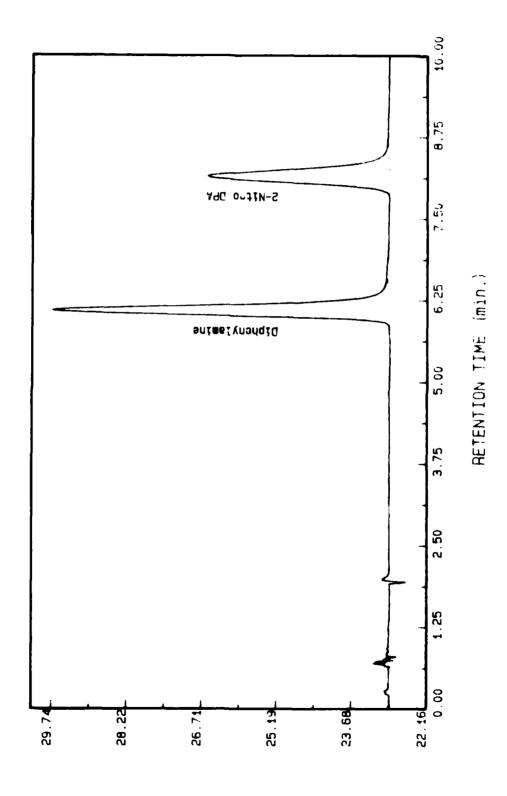


Figure 8. Typical responses of LC/EC detector at 1.05 V.

TABLE 3. HYDRA 70 2-NDPA Analysis by LC/FC 3 1.05 V

Results	SAMPLE !	SAMPLE 2
	2.0181 percent	2.0419 persent
	1.9964 percent	1.9899 parsent
	1.9788 percent	1.9939 per ent
AVG	1.9979 percent	2.0086 per est
SD	0.0197	0.0289
RSD	0.985	1.411

IV. CONCLUSIONS

The EC detector has been demonstrated to be a sensitive detector that is selective for propellant stabilizers. Analytical chemists concerned only wire the stabilizer levels of a propellant can select a solvent that will element stabilizers rapidly being only concerned with resolution of the stabilizer set its daughter products. This technique will be most useful for routine same alance procedures.

The internal standards used in this study were stabilizers themselves such as MNA, DPA, and 4-NDPA. As a follow-on study a variety of materials should be evaluated as internal standards. The optimum class of compounds would have a good UV response at 214 nm or at 254 nm and have a good response at above +0.9 V vs Ag/AgCl reference electrode. The class of compounds should be commercially available in high purity, have good stability in air, water, and most organic solvents, and have enough members to provide a wide scleation of retention times.

DISTRIBUTION

	No. Copies
ं अञ्चलकार भी स्टब्स	
.S. Army Materiel Command	
ATTN: AMCOA, Mr. S. Lorber	1
5 PM Sysenbower Avenue	
Nlexasteia, VA - 22333~0001	
Commander	
7.5. Army Materiel Command	
ATTN: AMCA, Mr. J. Inckerd	1
SPC] Fisenhower Avenue	
Alexandria, 74 (22333-000)	
To make the second seco	
.S. Brow Markers at Sommand	
ATTY AMORAT, Mr. F. 145ala	1
on the Company of the	
No. 2 x 4 x 4 x 5 x 4 x 5 x 5 x 3 x 3 x 3 x 4 x 6 x 6 x 6 x 6 x 6 x 6 x 6 x 6 x 6	
were growth and the second and the s	
。 、 April Materies () (mmail)	
1777 AM A-27, Mr. E. Shullo	•
Smith own to Avenue	
A service of the service of the Association of the service of the	
、 、 きがから (Made がこん) - 1 (MBBA) M	
ATTN: ATTN: ATTO ATTO MEN WE SERVED	1
Specifical Average	
10 - Control (1) - A - 22333 - 1903	
Logistation of the total of the community of the communi	
STORE AND AND AND ME. R. TEMPE	1
Control we convert Albert de	
The state of the s	
and the contract of the second	
STORE STATE OF STATE	1
A CONTRACTOR OF THE CONTRACTOR	
The state of the s	
er y Core	
At we March as well the following	
Research of the secretary secretary secretary	
ATTEL SMOUN-PRMATE, Mr. S. Karlin	1
1. Sept. 1. 1. 1. 1. Sept. 1.	

	No. Copies
Commander U.S. Army Materiel Command	
ATTN: AMCQA-ER, Mr. L. Garcia-Baco	1
5001 Eisenhower Avenue	•
Alexandria, VA 22333-0001	
The name of the same of the sa	
Commander	
U.S. Army Materiel Command	
ATTN: AMCTM-S, Mr. T. Wolf	1
5001 Eisenhower Avenue	
Alexandria, VA 22333-0001	
Comments	
Commander U.S. Army Materiel Command	
ATTN: AMCMT-M, Mr. G. Wagner	1
5001 Eisenhower Avenue	•
Alexandria, VA 22333-0001	
nichalidita, the 22775 oool	
Commander	
U.S. Army Armament, Munitions	
and Chemical Command	
ATTN: AMSMC-QA(R), Mr. L. Griffin	1
Rock Island, IL 61299-6000	
Commander	
".S. Army Armament, Munitions	
and Chemical Command	
ATTN: AMSMC-QAK-B(R), Mr. R. Fer	1
Rock Island, IL 61299-6000	
Commander	
S. Army Armament, Munitions	
and Chemical Command	!
ATTN: AMSMC+OA(D), Mr. H. Lazar Dover, NJ 07801-5001	•
Commander	
C.S. Army Armament, Munitions	
and Chemical Command	
ATTN: AMSMC-DAS-T(D), Mr. Glacobbe	:
Pover, NJ 07801-5001	
Commandad	
Commander	
and Chemical Command	
ATTN: SMCPB-OAL, Mr. V. Warren	
Pine Bluff Arsenal	
Pine Bluff, AR 71602-9500	
the same of exercises the contract of the cont	

	No. Copies
Commander	
U.S. Army Armament, Munitions	
and Chemical Command	
ATTN: AMSMC-SCM-P(D), Mr. J. Beetle	1
Dover, NJ 07801-5000	•
100ver, 45 0700r 7000	
Commander	
U.S. Army Jefferson Proving Ground	
ATTN: STEJP-TD, Mr. I. K. Peterson	1
Madison, IN 47250	-
Commander	
F.S. Army Yuma Proving Ground	
ATTN: STEYP-MTD, Mr. W. E. Brooks	1
Yama, 47 85365	
ਿਆmander	
7.5. Army White Sands Missile Range	
ATTN: STEWS-OA, Mr. C. Treat	1
White Sands Missile Range, VM 88002	
• •	
i ommander	
1.S. Army Tropts Test Ctr.	
ATTN: STRICHLINEM	1
(\mathbf{A},\mathbf{P}) .	
Mi 171 - 14(H)4	
·	
of ommander	
S. Army Troop Support & Aviation	
Materiel Readiness Command	
ATIN: AMSTS-O, Mr. W. G. Creel	1
•3(h) Goodfellow Boulevard	
St. 'outs, MO 63120-5720	
€ nomander	
Anniston Army Depot	
ATTN: SBSAN-DDA, Mr. Pennington	1
Anniston, Al. 36201	
·	
√ommanter	
C.S. Army Depot System Command	
ATTN: AMSDS-OM, Mr. R. Newman	1
Thumbersbury, PA 17201-4170	
√ommanter	
Corpus Christi Army Depot	
ATTN: SDSCC-O, Mr. L. F. Frvin	1
Corpus Christi, TV 78419	

THE RESERVE OF THE PARTY OF THE

	No. Copies
Commander Letterkenny Army Depot ATTN: SDSLE-Q, Mr. W. J. Walker Chambersburg, PA 17201-4170	1
Commander New Cumberland Army Depot ATTN: SDSNC-Q, Mr. A. T. Holderbach New Cumberland, PA 17070	1
Commander Red River Army Depot ATTN: SDSRR-Q, Mr. W. D. Wuertz Texarkana, TX 75501	1
Commander Lexington-Blue Grass Activity ATTN: SDSAN-LQ, Mr. P. Jackson Lexington, KY 40511	1
Commander Sacramento Army Depot ATTN: SDSSA-Q, Mr. C. H. Clawson Sacramento, CA 95801	1
Commander Savanna Army Depot ATTN: SDSSV-VE, LTC P. E. James Savanna, IL 61074	1
Commander Seneca Army Depot ATTN: SDSSE-R, Mr. P. W. Chavez Romulus, NY 14541	1
Commander Sharpe Army Depot ATTN: SDSSH-Q, Mr. J. E. Seyfried Lathrop, CA 95330	1
Commander Sierra Army Depot ATTN: SDSSI-QA, Mr. L. B. Lampe Herlong, CA 96113	1
Commander Tooele Army Depot ATTN: SDSTE-QA, Mr. R. M. Rich Tooele UT 84074	1

	No. Copies
Commander	
Pueble Depot Activity	
ATTN: SDSTE-PUG, Mr. P. S. Coen	1
Pueblo, CO 81001	
Commander	
Tobyhanna Army Depot	•
ATTN: SDSTO-Q, Mr. W. J. Lord	1
Tobyhanna, PA 18466	
Director	
U.S. Army Industrial Base	
Engineering Activity	
ATTN: AMXIB-MT, Mr. Del Brim	1
Rock Island, IL 61299-7260	-
Noch Island, In Olly, 1800	
Commander	
Watervliet Arsenal	
ATTN: SMCWV-QA, Mr. J. Miller	1
Watervliet, NY 12189-5000	
Commander	
Watervliet Arsenal	_
ATTN: SMCWV-QAE-S, Mr. S. Krupski	1
Watervliet, NY 12189-5000	
Commander	
Watervliet Arsenal	
ATTN: SMCWV-PPI, Mr. P. Phelan	1
Watervliet, NY 12189-5000	•
waterviret, Nr. 12107 3000	
Commander	
U.S. Army Armament, Munitions	
and Chemical Command	
ATTN: AMSMC-QAS-A(D), Mr. Geza Pap	1
Dover, NJ 07801-5000	
Commander	
U.S. Army Armament, Munitions	
and Chemical Command	,
ATTN: AMSMC-QAC-E(E), Dr. W. J. Maurits	1
Aberdeen Proving Ground, MD 21010-5423	
Commander	
II.S. Army Natick R&D Center	
ATTN: STRNC-EML, Mr. R. Merullo	1
Kansas Street	
Nortick MA 01760-5014	

	No. Copies
Commander	
U.S. Army Tank-Automotive Command ATTN: AMSTA-Q, Mr. L. Barnett Warren, MI 48090-2004	1
Commander	
U.S. Army Tank-Autmotive Command ATTN: AMSTA-QAT, Mr. D. Gamache Warren, MI 48090-2004	-
Commander	
U.S. Army Tank-Automotive Command ATTN: AMSTA-QAT, Mr. F. Braun Warren, MI 48090-2004	1
Commander	
U.S. Army Tank-Automotive Command ATTN: AMSTA-RCK, Dr. J. Chevalier Warren, MI 48090-2004	1
Commander	
U.S. Army Tank-Automotive Command ATTN: AMSTA-RCKM, Mr. D. Ostberg Warren, MI 48090-2004	1
Commander	
U.S. Army Test & Evaluation Command	_
ATTN: AMSTE-AD-M, R. G. Sheldon Aberdeen Proving Ground, MD 21005	1
Commander U.S. Army Test & Evaluation Command	
ATTN: AMSTE-AD-M, Mr. R. Huddleston	1
Aberdeen Proving Ground, MD 21005	
Commander	
U.S. Army Test & Evaluation Command ATTN: AMSTE-AD-M, Mr. W. Deaver	1
Aberdeen Proving Ground, MD 21005	
Commander	
U.S. Army Test & Evaluation Command	1
ATTN: AMSTE-AD-M, Dr. N. Pentz Aberdeen Proving Ground, MD 21005	1
Commander	
Aberdeen Proving Ground	
ATTN: AMSTE-MT-G, Mr. C. Klarich	1
Aberdeen Proving Ground, MD 21005	

	No. Copies
Commander U.S. Army Electronic Proving Ground ATTN: STEE-PMT, COL G. L. Kosmider Fort Huachuca, AZ 85613	1
Commander U.S. Army Dugway Proving Ground ATTN: STEDP-P0-Q Dugway, UT 84022	1
Commander U.S. Army Aviation System Command ATTN: AMSAV-Q, Mr. E. J. Hollman St. Louis, MO 63120-1798	1
Commander U.S. Army Aviation System Command ATTN: AMSAV-QE, Mr. J. T. Conroy St. Louis, MO 63120-1798	1
Commander U.S. Army Aviation System Command ATTN: AMSAV-QE, Mr. A. Spratt St. Louis, MO 63120-1798	1
Commander U.S. Army Depot Systems Command ATTN: AMSDE-Q, Mr. W. D. Henry Letterkenny Army Depot Chambersburg, PA 17201-4170	1
Commander U.S. Army Communications-Electronics Command ATTN: AMSEL-PA, Mr. A. D'Angelo Fort Monmouth, NJ 07703-5023	1
Commander ".S. Army Communications=Electronics Command ATTN: AMSEL=PA=M, Mr. F. Di Giovanni Fort Mondouth, NJ 07703-5023	1
Commander U.S. Army Electronics R&D Command ATTN: AMDEL-PA, Mr. S. Alster 2800 Powder Mill Road Adelphi MD 20783-1197	1

	No. Copies
Commander	
U.S. Army Electronics R&D Command	
ATTN: AMDEL-PA, Mr. J. Goon	1
2800 Powder Mill Road	
Adelphi, MD 20783-1197	
Commander	
U.S. Army Night Vision &	
Electro-Optics Laboratory	
ATTN: DELNV-SE, R. Stefanik	1
Fort Belvoir, VA 22060	
Commander	
U.S. Army Belvoir R&D Center	
ATTN: STRBE-HE, Ms. V. Lee	1
Fort Belvoir, VA 22060-5606	
Commander	
U.S. Army Belvoir R&D Center	
ATTN: STRBE-HE, Mr. K. Harris	1
Fort Belvoir, VA 22060-5606	
Commander	
U.S. Army Armament, Munitions	
and Chemical Command	
ATTN: AMSMC-QAR-I-(D),	
Mr. M. H. Weinberg	1
Dover, NJ 07801-5000	
Director	
Army Materials and Mechanics	
Research Center	
ATTN: AMXMR-D	1
AMXMR−DD AMXMR−M	1
AMXMR -0	1 1
AMXMR-OP, G. Hagnauer	1
AMXMR-P	i
AMXMR-S	1
AMXMR-ST	ī
AMXMR-STN	ī
Watertown, MA 02172-0001	-
U.S. Army Materiel Systems Analysis	
Activity	
ATTN: AMXSY-MP	1
Abardson Proving Cround MD 21005-5071	

	No. Copies
ITT Research Institute	
ATTN: GACIAC	1
10 W. 35th Street	
Chicago, IL 60616	
AMSMI-RD, Dr. McCorkle	1
Dr. Rhoades	1
AMSMI-RD-PR-M, Mr. J. Wright	1
-T, Mr. B. J. Alley	1
AMSMI-RD-CS-T	1
AMSMI-RD-CS-R	15
AMSMI-GC-IP, Mr. Bush	1
AMSMI-SE-MT, Mr. J. Jackson	1
AMSMI-SE-MT, Mr. L. Chapman	1
AMSMI-QA-QE, Mr. T. McVey	1
-QA, Mr. T. Howard	1

CONTRACTOR OF STREET